

Standardised and reproducible analysis of mass spectrometry-based single-cell proteomics data

Replication of the SCoPE2 analysis by Specht et al. 2019

Christophe Vanderaa, Laurent Gatto

Computational Biology Unit (CBIO)
de Duve Institute
UCLouvain
Belgium

18 August 2020

Outline


► Expectation



► Expectation



► Reality

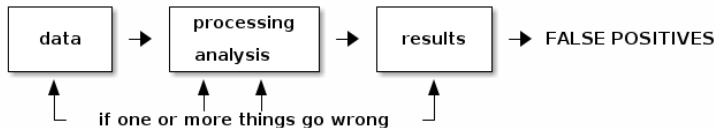

- ▶ **Replication-based development** agreement between the developer, the data producer and the user.

Replication is the first step to define sound data infrastructure and principled analysis.

- ▶ SCoPE2 quantifies thousands of proteins x thousands single-cells
- ▶ Full protocole available
- ▶ Full analysis script and data available

- ▶ Contribute a standardised and principled data and analysis that is broadly applicable.
- ▶ Open, transparent and reproducible computational infrastructure to further improve data analysis and interpretation.
- ▶ R and Bioconductor (?) offer an ideal environment to attain these goals.


Implemented in the `scp` package.

Outline

Data infrastructure: QFeatures¹

scp package

QFeatures: data framework dedicated to manipulate and process MS-based quantitative data.


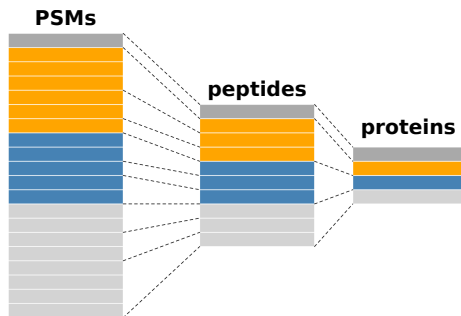


¹Gatto (2020)

Data infrastructure: QFeatures¹

scp package

QFeatures: data framework dedicated to manipulate and process MS-based quantitative data.

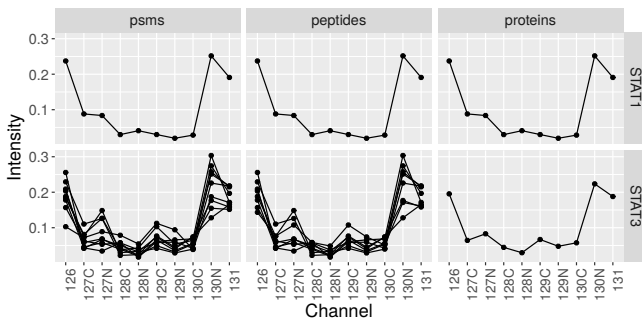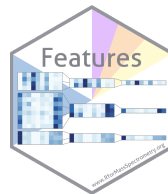


¹Gatto (2020)

Data infrastructure: SingleCellExperiment^{2,3}

scp package



`SingleCellExperiment`: provides dedicated framework for single-cell data analysis.

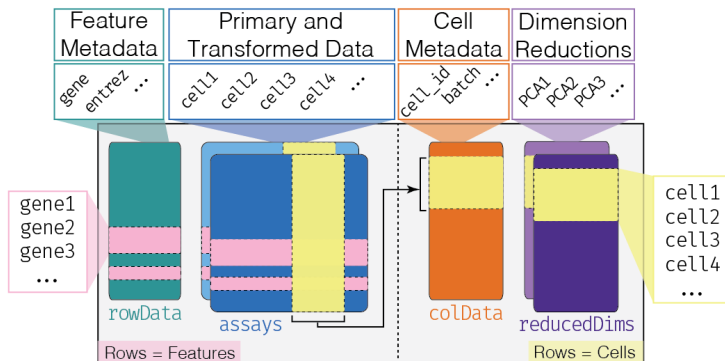


SingleCellExperiment

²Lun and Risso (2020)

³Amezquita et al. (2019)

`scp` = `SingleCellExperiment` + `QFeatures`

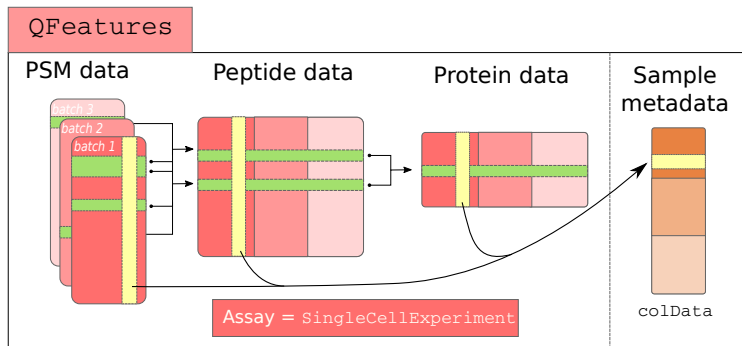

AND functions dedicated to processing and analyzing SCP data.

Load the SCoPE2 dataset called `specht2019v2`⁴

```
1 library(scpdata)
2 data("specht2019v2")
```

Dataset overview

```
1 show(specht2019v2)
```

```
An instance of class QFeatures containing 179 assays:
[1] 190222S_LCA9_X_FP94AA: SingleCellExperiment with 2823 rows and 11 col...
[2] 190222S_LCA9_X_FP94AB: SingleCellExperiment with 4297 rows and 11 col...
[3] 190222S_LCA9_X_FP94AC: SingleCellExperiment with 4956 rows and 11 col...
...
[177] 191110S_LCB7_X_APNOV16plex2_Set_9: SingleCellExperiment with 4626 r...
[178] peptides: SingleCellExperiment with 9208 rows and 1018 columns
[179] proteins: SingleCellExperiment with 2772 rows and 1018 columns
```

Tabular data (generated by MaxQuant, ProteomeDiscoverer, ...) can be converted to `QFeatures` using the `readSCP()` function.

⁴Specht et al. (2019)

Sample metadata common to **all assays** are stored in a single table, the `colData`

Set	Channel	SampleType	lcbatch	sortday	digest
190222S_LCA9_X_FP94AA	RI1	Carrier	LCA9	s8	N
190222S_LCA9_X_FP94AA	RI2	Reference	LCA9	s8	N
190222S_LCA9_X_FP94AA	RI3	Unused	LCA9	s8	N
190222S_LCA9_X_FP94AA	RI4	Macrophage	LCA9	s8	N
190222S_LCA9_X_FP94AA	RI5	Macrophage	LCA9	s8	N
190222S_LCA9_X_FP94AA	RI6	Macrophage	LCA9	s8	N
...

1. Load data

PSM data

```
[1] 190222S_LCA9_X_FP94AA: SingleCellExperiment with 2823 rows and 11 columns
[2] 190222S_LCA9_X_FP94AB: SingleCellExperiment with 4297 rows and 11 columns
[3] 190222S_LCA9_X_FP94AC: SingleCellExperiment with 4956 rows and 11 columns
...
[177] 191110S_LCB7_X_APNOV16plex2_Set_9: SingleCellExperiment with 4626 rows and 16 columns
```

1. Load data

PSM data

```
[1] 190222S_LCA9_X_FP94AA: SingleCellExperiment with 2823 rows and 11 columns
[2] 190222S_LCA9_X_FP94AB: SingleCellExperiment with 4297 rows and 11 columns
[3] 190222S_LCA9_X_FP94AC: SingleCellExperiment with 4956 rows and 11 columns
...
[177] 191110S_LCB7_X_APNOV16plex2_Set_9: SingleCellExperiment with 4626 rows and 16 columns
```

2. PSM filtering
3. Expression channel by reference channel division
4. PSM to peptides aggregation
5. Joining sets
6. Single cells filtering based on median CV
7. Normalization
8. Removal of highly missing peptides
9. Log-transformation

Peptide data

```
[178] peptides: SingleCellExperiment with 9208 rows and 1018 columns
```


1. Load data

PSM data

```
[1] 190222S_LCA9_X_FP94AA: SingleCellExperiment with 2823 rows and 11 columns  
[2] 190222S_LCA9_X_FP94AB: SingleCellExperiment with 4297 rows and 11 columns  
[3] 190222S_LCA9_X_FP94AC: SingleCellExperiment with 4956 rows and 11 columns  
...  
[177] 191110S_LCB7_X_APNOV16plex2_Set_9: SingleCellExperiment with 4626 rows and 16 columns
```

2. PSM filtering
3. Expression channel by reference channel division
4. PSM to peptides aggregation
5. Joining sets
6. Single cells filtering based on median CV
7. Normalization
8. Removal of highly missing peptides
9. Log-transformation

Peptide data

```
[178] peptides: SingleCellExperiment with 9208 rows and 1018 columns
```

9. Peptides to proteins aggregation
10. Normalization
11. Imputation
12. Batch correction

Protein data

```
[179] proteins: SingleCellExperiment with 2772 rows and 1018 columns
```

1. Load data

PSM data

```
[1] 190222S_LCA9_X_FP94AA: SingleCellExperiment with 2823 rows and 11 columns  
[2] 190222S_LCA9_X_FP94AB: SingleCellExperiment with 4297 rows and 11 columns  
[3] 190222S_LCA9_X_FP94AC: SingleCellExperiment with 4956 rows and 11 columns  
...  
[177] 191110S_LCB7_X_APNOV16plex2_Set_9: SingleCellExperiment with 4626 rows and 16 columns
```

2. PSM filtering

3. Expression channel by reference channel division

4. PSM to peptides aggregation

5. Joining sets

6. Single cells filtering based on median CV

7. Normalization

8. Removal of highly missing peptides

9. Log-transformation

Peptide data

```
[178] peptides: SingleCellExperiment with 9208 rows and 1018 columns
```

9. Peptides to proteins aggregation

10. Normalization

11. Imputation

12. Batch correction

Protein data

```
[179] proteins: SingleCellExperiment with 2772 rows and 1018 columns
```

Outline

Filter out features based on the feature metadata

Example: filter out reverse hits. The filter is applied to the `Reverse` field in the feature metadata

```
1 filterFeatures(specht2019v2 ,  
2               ~ Reverse != "+")
```

Source code in `QFeatures`

Some QC metrics are not computed by MaxQuant:

- ▶ Sample to carrier ratio: discard samples with intensities higher than expected
- ▶ Peptide FDR: expected proportion of features wrongly assigned to a given peptide
- ▶ Cell median coefficient of variation: reliability of the protein quantification in a cell.


Example:

```
1 computeMedianCV(specht2019v2,
2                 i = "peptides",
3                 proteinCol = "protein",
4                 peptideCol = "peptide",
5                 batchCol = "Set")
```

Source code in scp


QC metrics are stored in the data set for plotting or subsetting

```
1 library(tidyverse)
2 specht2019v2[["peptides"]] %>%
3   colData %>%
4   data.frame %>%
5   ggplot(aes(x = MedianCV,
6             fill = SampleType)) +
7   geom_histogram() +
8   geom_vline(xintercept = 0.4)
```

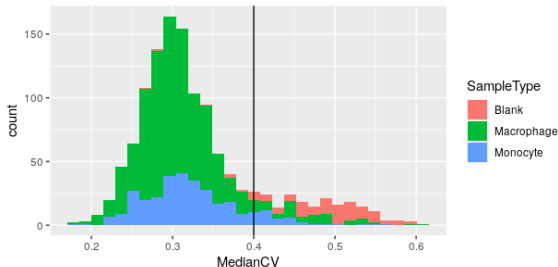

Feature aggregation = combine features into a higher-level structure.

Example: aggregate peptides to proteins

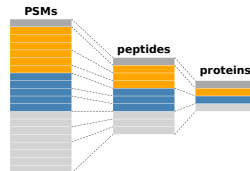

- ▶ Combine the quantitative data from multiple peptides to a single protein
- ▶ Store the relationship between the protein and the aggregated peptides

```
1 aggregateFeatures(specht2019v2 ,  
2                   i = "peptides",  
3                   name = "proteins",  
4                   fcol = "protein",  
5                   fun = colMedians, na.rm = TRUE)
```

Source code in `QFeatures`

0's can be either **biological** or **technical** zero. They are better related by NA's.

```
1 zeroIsNA(specht2019v2 ,  
2           i = "peptides")
```

Remove highly-missing features (e.g. $\geq 99\%$)

```
1 filterNA(specht2019v2 ,  
2           i = "peptides" ,  
3           pNA = 0.99)
```

Impute missing data

```
1 impute(specht2019v2 ,  
2         i = "proteins" ,  
3         method = "knn" ,  
4         k = 3)
```

Source code in `QFeatures`

Common data transformation can easily be applied such as **log-transformation** or **normalization**.

Example: \log_2 -transformation:

```
1 logTransform(specht2019v2,  
2             i = "peptides",  
3             base = 2,  
4             name = "peptides_log")
```

Source code in `QFeatures`

Custom function can be applied to the data set, for example batch correction using 'ComBat'. Three-step procedure:

1. Extract the assay to process

```
1 x <- specht2019v2[["proteins"]]
```



2. Apply the custom function

```
2 assay(x) <- ComBat(assay(x), ...)
```


3. Add the processed assay in the dataset

```
3 addAssay(specht2019v2, x, name = "proteins_batch_corrected")
```


Outline

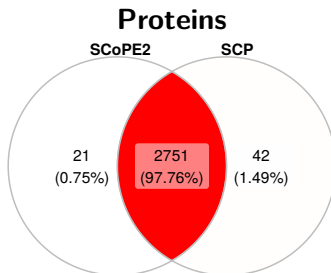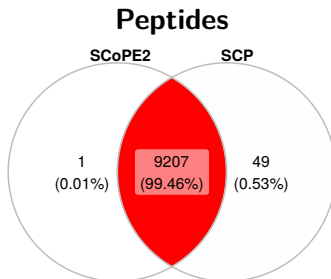

Peptides

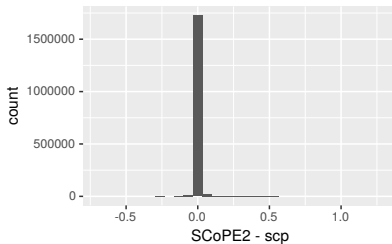



Proteins

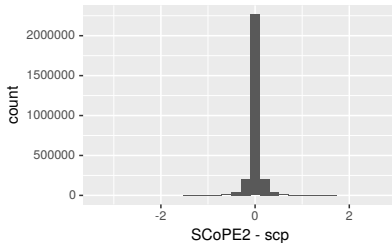

Weighted PCA on the **protein** data

SCoPE2 (Figure 3a in preprint)

scp

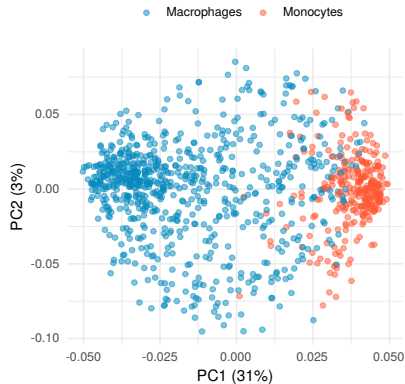
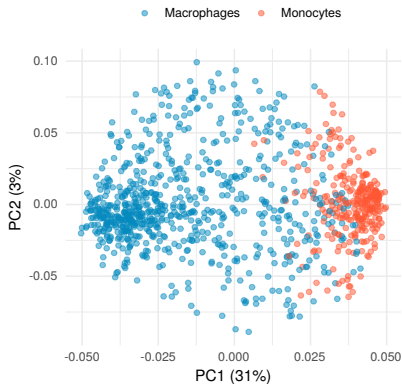


Note about reproducibility

Replication results

Good software development includes continuous maintenance and improvement **BUT** might impact reproducibility



Example: batch correction using 2 versions of the **ComBat** algorithm
(**sva** package)

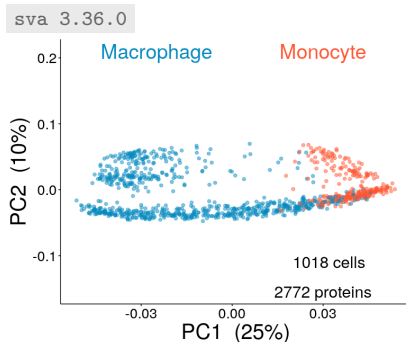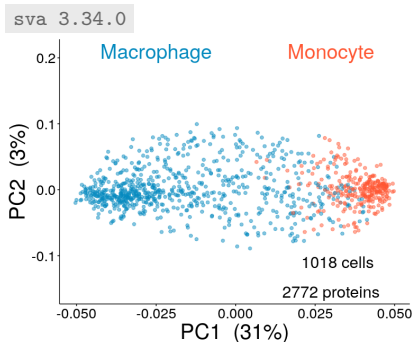

Documenting software version is **essential** for reproducible work

Outline

- ▶ MS-based single cell proteomics: young field, with many challenges and great progress. `scp` to address the need for principled and reproducible data analysis.
- ▶ `scp` isn't specific to SCoPE2/TMT data, applicable to other LF protocols such as nanoPOTS (Williams et al. (2020); Cong et al. (2020)).
- ▶ `scp` and `SingleCellExperiment`: same infrastructure for single cell proteomics and RNA sequencing.
- ▶ Tool for novel computational developments.

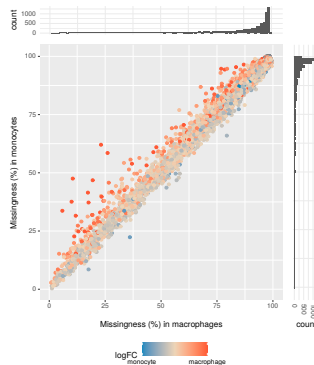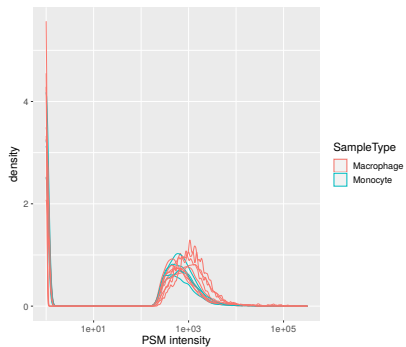

Resources

- ▶ `scp`: <http://UClouvain-CBIO.github.io/scp>
- ▶ `scpdata`: coming soon
- ▶ `QFeatures`: <http://rformassspectrometry.org>
- ▶ `SingleCellExperiment`: Bioconductor
- ▶ **Slides**: <http://bit.ly/2020SCP> under CC-BY SA

Resources

- ▶ `scp`: <http://UClouvain-CBIO.github.io/scp>
- ▶ `scpdata`: coming soon
- ▶ `QFeatures`: <http://rformassspectrometry.org>
- ▶ `SingleCellExperiment`: Bioconductor
- ▶ **Slides**: <http://bit.ly/2020SCP> under CC-BY SA

Acknowledgements

- ▶ Nikolai Slavov, Harrison Specht, Ed Emmott.
- ▶ Fonds National de la Recherche Scientifique (FNRS)
- ▶ **Thank you for your attention**

References I